

Portable and durable.  
For visiting nurses, pediatric  
offices and family practices.

Lightweight | Durable  
Highly accurate | Battery and AC

Doran. Powering Performance.

doran®

800.264.4107 | [doranmedical.com](http://doranmedical.com)

**DS4500 INFANT  
LACTATION SCALE**

45lb (20kg) capacity with  
weighing sensitivity to  
track even the smallest  
weight changes. Advanced  
best-in-class features  
include one-button Feed  
function for recording  
before and after feeding  
and amount fed.  
Stylish design with  
one-piece construction  
for easy sanitation.  
Affordably priced.



# Infant Lactation Scale


The DS4500 is part of a complete line of high-quality Doran medical scales, designed to be durable and highly accurate, with innovative safety and performance features that support patient comfort and easy operation. Backed by Doran three-year warranty.

## INNOVATIVE DESIGN AND PERFORMANCE FEATURES

- Weight is quickly locked-on and displayed with Doran exclusive Weight Lock technology.
- Large platform measuring 25" L x 11" D with curved sides for baby safety.
- Lightweight and portable with exceptional build quality.
- Battery or AC power for continuous operation.
- Audio signals indicate scale status and FEED operation status, allowing attention to be focused on the patient, rather than the weight display.
- Reweigh feature allows the infant to be weighed again without removing the baby from the scale.
- Quickly switch between lb, lb:oz or kg.
- EMR ready with USB port.
- Rubber feet to maintain stability.

## SPECIFICATIONS

<b>Capacity:</b>	45 lb x 0.1 oz, 45 x 0.005 lb, 20 kg x 2 g
<b>Accuracy:</b>	+/- 0.2 oz, 0.01 lb, 5 g
<b>Dimensions:</b>	25" x 13.75" W x 3.75" H
<b>Power:</b>	AC transformer included. 4 AA batteries provides up to one-year of weighments
<b>Display:</b>	High visibility 0.66" LCD display
<b>Measuring tape:</b>	Included. May be attached to the scale where it will be most convenient.
<b>Option:</b>	Carrying case

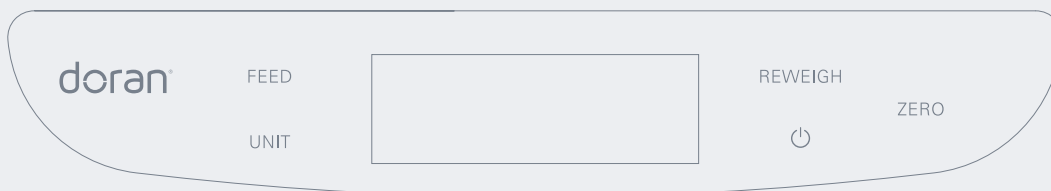

## FIVE-BUTTON CONTROLS

### Power

On/Off with adjustable Automatic Shutoff Timer for battery power applications. It is disabled when in feed mode.

### Feed

For taking weigh measurements before and after infant feeding and determining amount fed.

### Reweigh

Button to reinitiate the weighing process without the patient needing to get off the scale platform.

### Unit

Any unit can be locked out to ensure accurate weighments.

### Zero

Feature can account for clothing or blankets.

